Action Sports Video App

RMG Enterprises, Inc

DBA Sportsvet, Inc (in South Carolina)

Personal Activity and Sports Analysis System

Action Analysis for Anyone


- Easily Capture Video
- Playback Video
 - Controlled Playback
 - Slow Motion
 - Stop or Pause Playback
- Review or Analyze Action Video
- Present to Athlete, Patient or Client
- Prepare a Report

Personal Analysis System

- An App that can be easily available and usable to capture onsite actions, activities or sports
- Phone or Tablet
 - Tablets have a larger viewing screen
 - Tablets allow for easy presentation to others
 - Phones are more readily available
- Affordable

Membership screen

- On this screen, you can purchase a membership plan.
- You also get a 7-day free trial.
- This subscription provides several features, including changing playback speed, editing snapshots, trimming, downloading, and applying filters to videos.
- This screen includes links to the Privacy Policy and Terms of Service, allowing you to visit them anytime.
- It also has a **Rate Us** option, so you can directly rate the app on the **Play Store** or **App Store**.

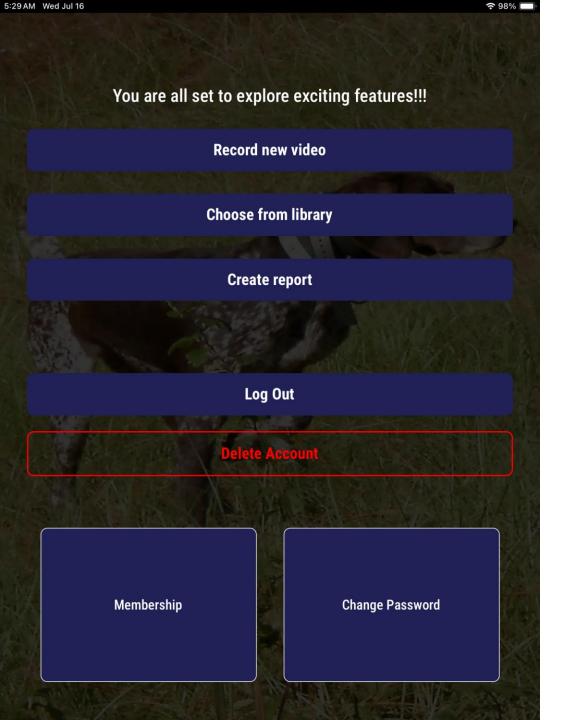


Using the Action Sports Video App


At the start you have three options


- Record new video
 - Record a new video for analysis
- Chose from library
 - Chose a previously recorded video for analysis
- Create report
 - Create a report about the video analysis for a patient or client

Using the Action Sports Video App

If you are viewing in "Landscape" mode, you may need to scroll down to see options at the bottom of window.

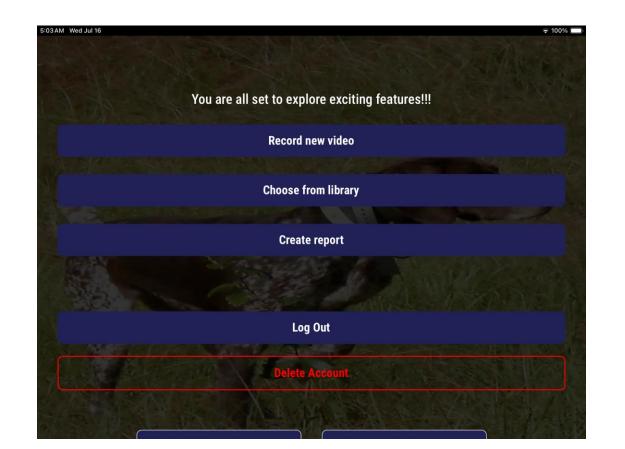


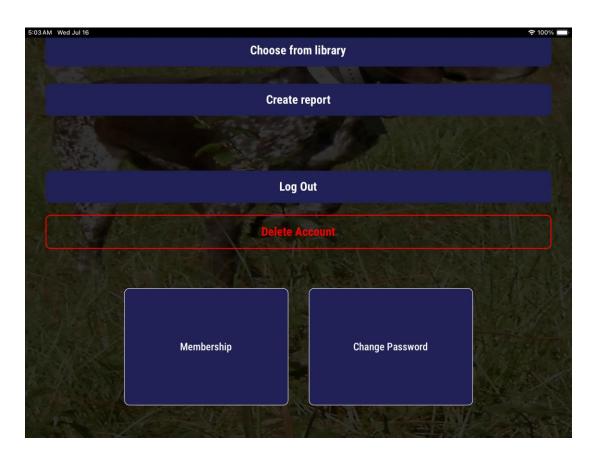

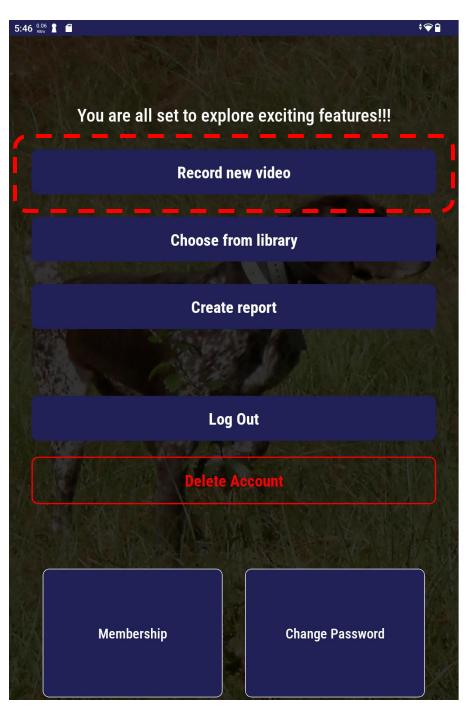

Record New Video


To record a new video, select "record new video"

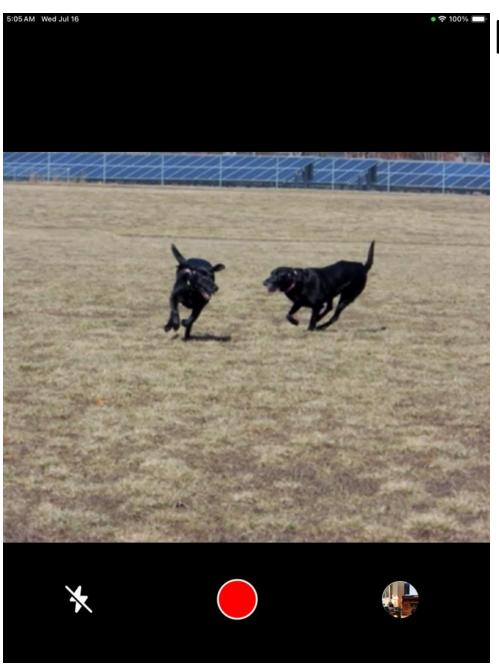

Record New Video

To record a new video, select "record new video"

- You can record the video in "Portrait" Mode
 - Many "Human" activities can be capture in this mode
- You can record in "Landscape" Mode
 - Many "Animal" activities can be captured in this mode

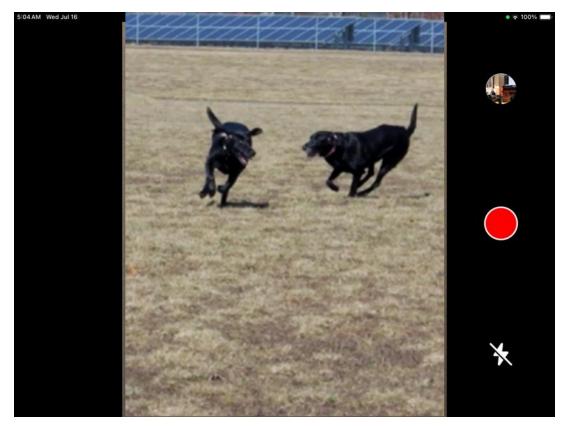

← Video Play

Delete

Edit



Video Play Screen


After the video is captured, the user is navigated to the video player screen.

At the top of this screen, there are three icons:

- **Delete**: Deletes the selected video.
- Edit: Navigates to the edit video screen.
- **Save**: Saves the video to the local device.


The video player includes the following functionalities:

- Rewind and forward by one frame at a time.
- Pause and resume playback.
- Extreme left button: Jumps to the start of the video.
- Extreme right button: Jumps to the end of the video
- You can also use the "slider" to move through the video

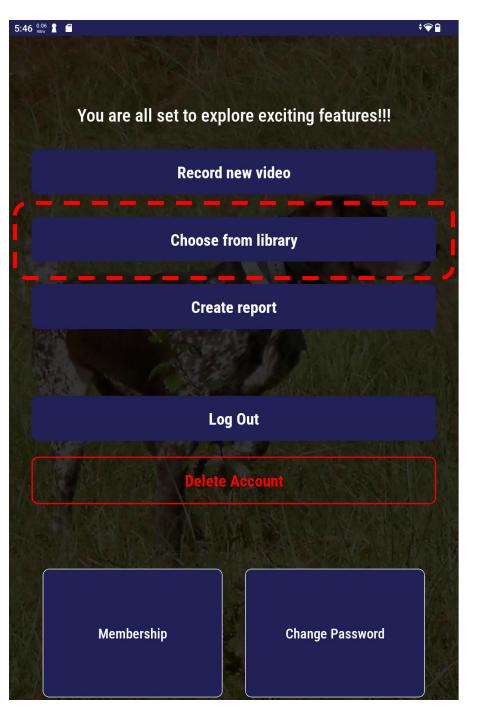

Chose from Library

To select a previously recorded video, select "Chose from Library"

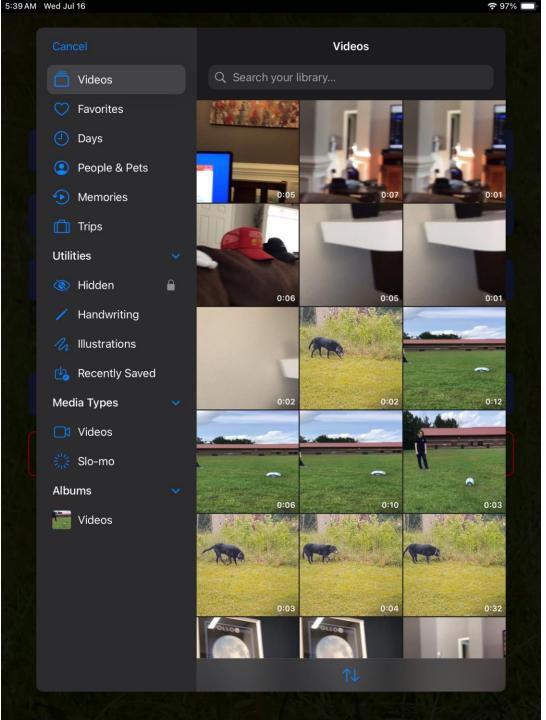

Selecting video from gallery

You can select any video from the library



Video Play Screen

After you select a video, this screen will appear.


The video player allows you to:

- Rewind and forward by one frame at a time.
- Pause and resume playback.
- Use the **extreme left** and **right** buttons to jump to the start and end of the video, respectively."
- You can also use the "slider" to move through the video
- Select the edit button to move to the next options




Editing Video Screen

- You can select any of the four options: Trim, Change
 Speed, Apply Filters, or Snapshot.
- You can either cancel the current video editing session or save the edited video by tapping the Save button.

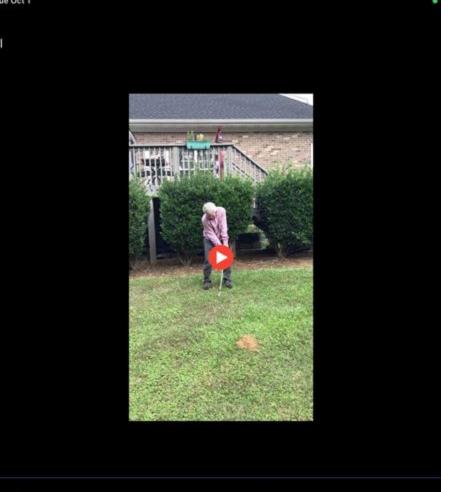



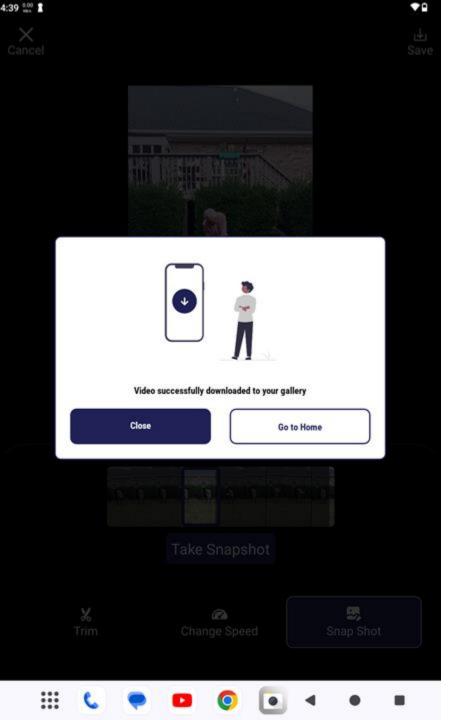


Video Successfully Download Popup

- After you have selected to save the video
- This pop-up confirms that your video has been successfully downloaded to your device.
- You can either close it and continue editing the video or go to the home screen.


Trim the video screen

- Here, you can trim the video by adjusting the thumbs from both sides and apply the trim by tapping the tick icon.
- You can either cancel the current video editing session or save the edited video by tapping the **Save** button.

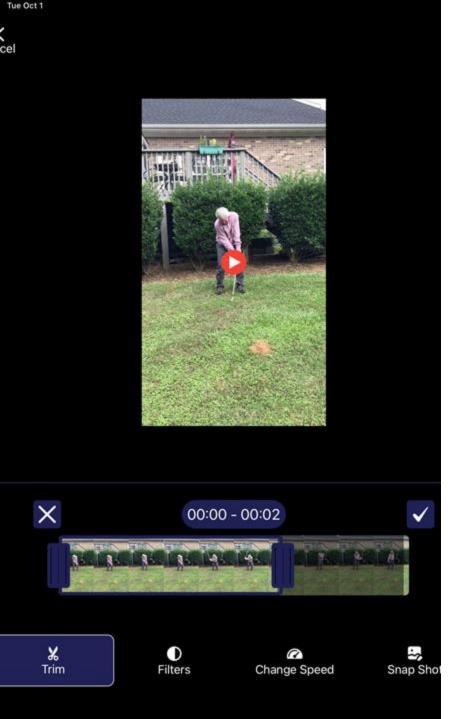

Trim Video Popup

- This pop-up confirms that your video has been trimmed.
- You can close the pop-up to continue editing the video.

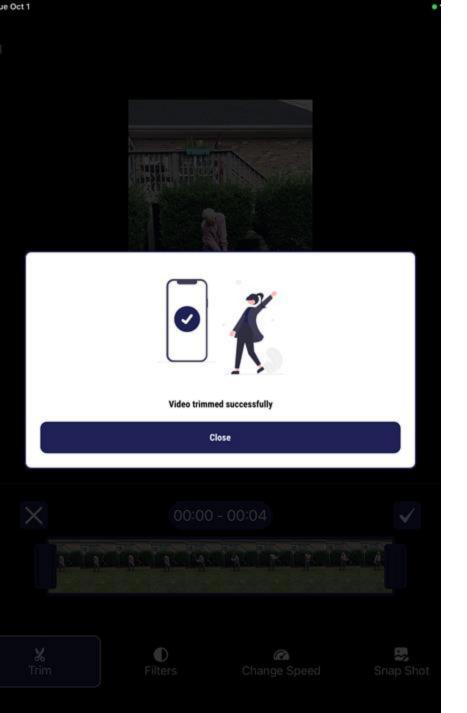

Apply Filters to the video screen

- You can select and apply any filter to the video.
- You can either cancel the current video editing session or save the edited video by tapping the Save button.







Tue Oct 1










Change Speed


Change Video Playback Speed

- Here, you can change the playback speed of the video being edited.
- After adjusting the speed, you can save the modified video to your device by tapping the Save button

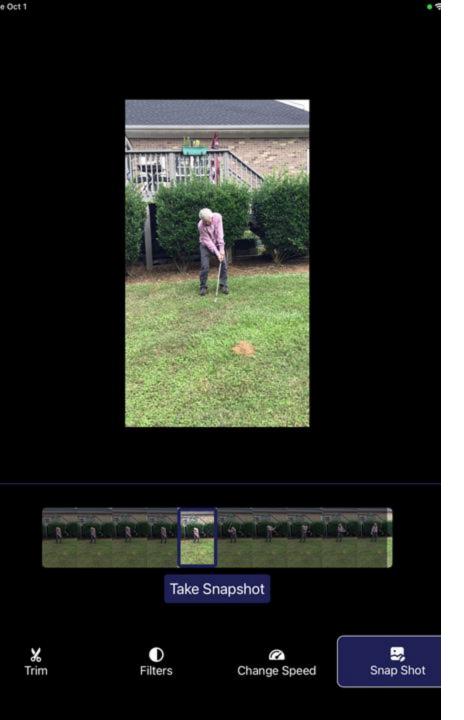


Snapshot of particular frame of video


- Here, you can select any frame by sliding through the video frames
- After selecting a frame, you can edit it on the next screen

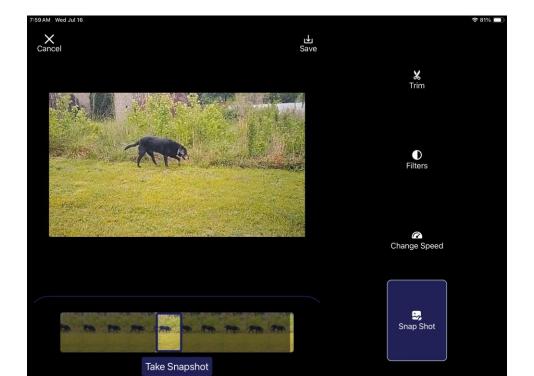


- Here, you will see the snapshot you have captured.
- If editing options are not visible, tap on the snapshot to see the editing options
- Note: Pen icon can be either at top or bottom of the screen based on device model.

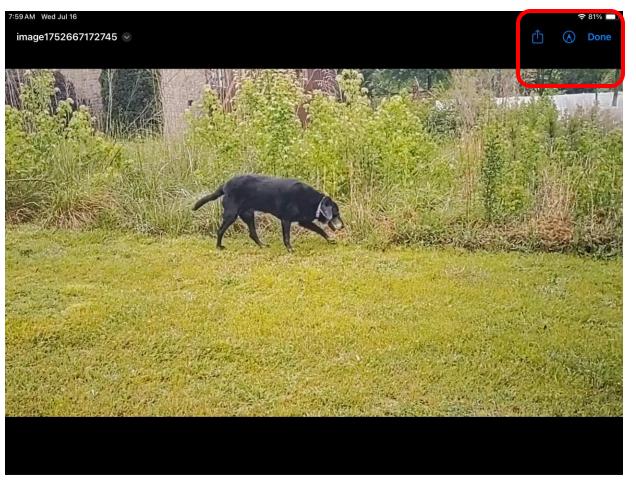




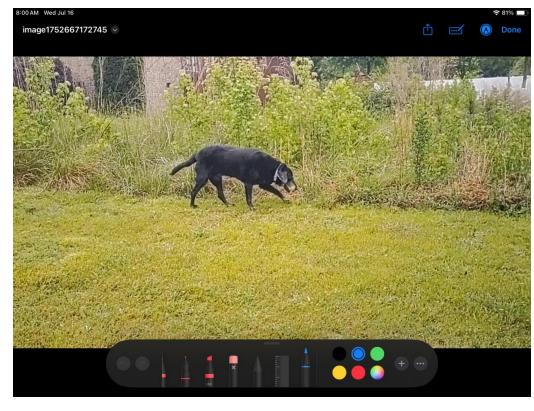



- Select any of the drawing utensils and draw shapes on the image
- You can select any color from the palette

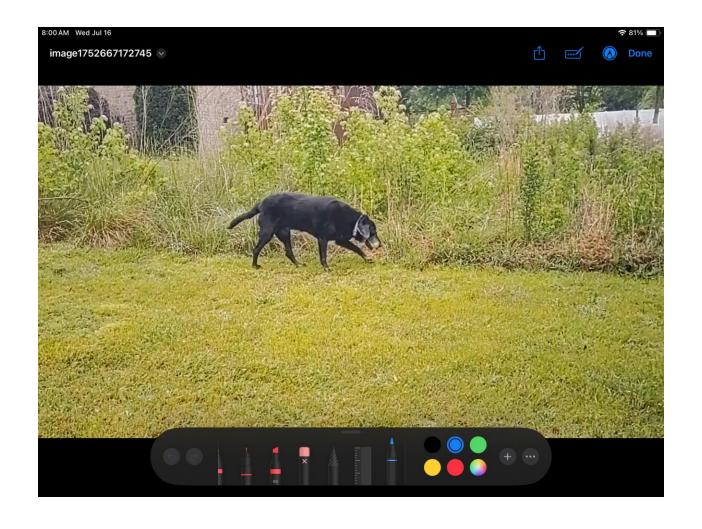


- To add shapes, text and other otions, tap on the "+" sign
- To add a shape, tap on 'Add Shape', select the shape and adjust its color.

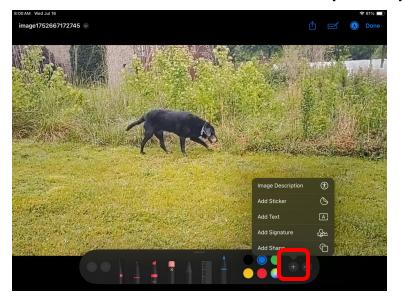

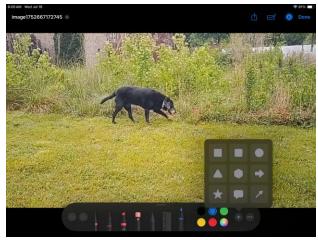

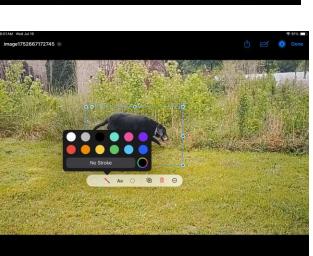





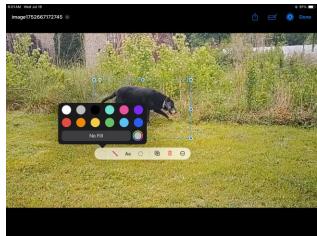



- To add Text and other options
- Here, tap on the + sign to open those options Tap on 'Add Text' and write your text

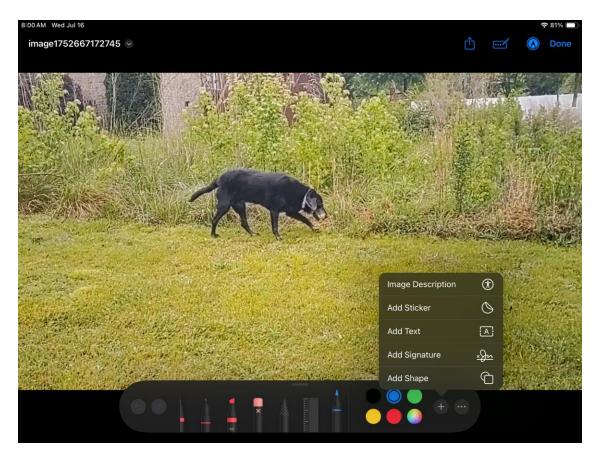


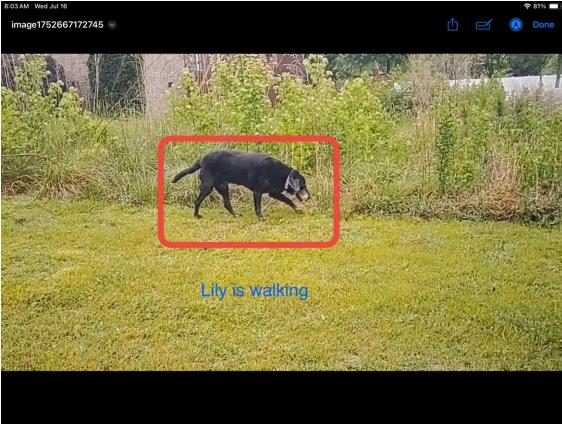

- Save the edited Snapshot by clicking the "Save" icon
- Select "Save" the edited snapshot into your device
 - This will save the image at a location the report can find it
- Note: Save option position can vary in different models of device
- Once you are done with editing, tap on 'Done'
- If you have already saved the image just select discard edits


Save Picture Options




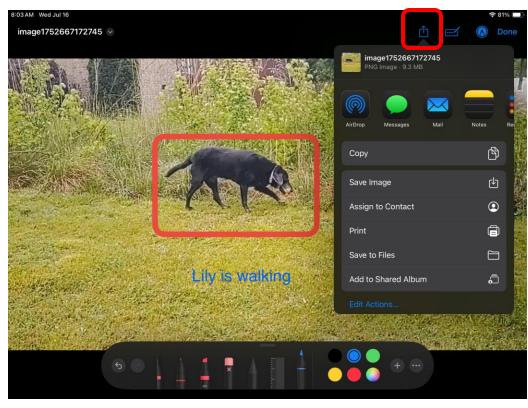

Snapshot successfully downloaded pop up


- When you tap the Save button, a pop-up will appear confirming that the snapshot has been successfully downloaded to your gallery.
- You can continue editing by tapping the Close button.

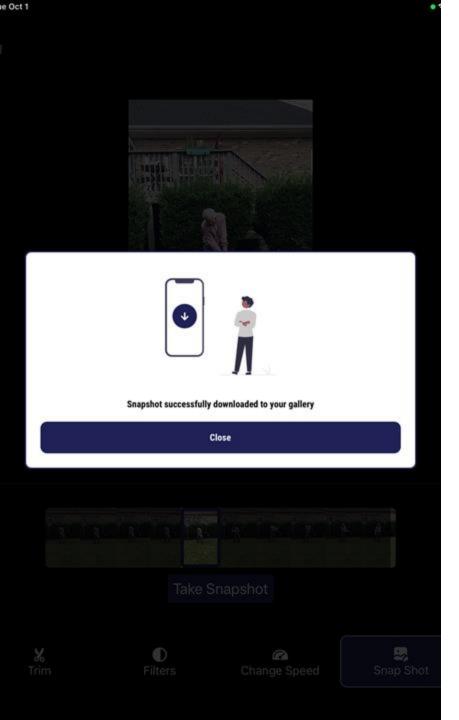

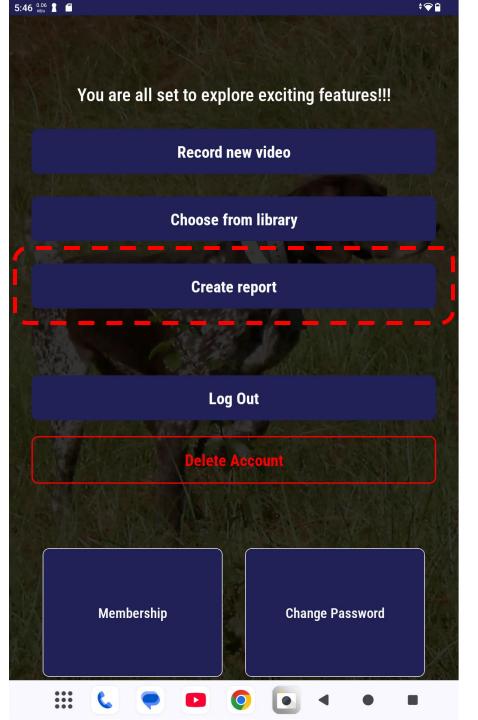

A report can be created to provide information to clients, customers or others who may have interest in the actions captured by the video.

- Initially, a blank screen will appear if no previous reports have been created
- Click the "Plus" button to begin creating a report.

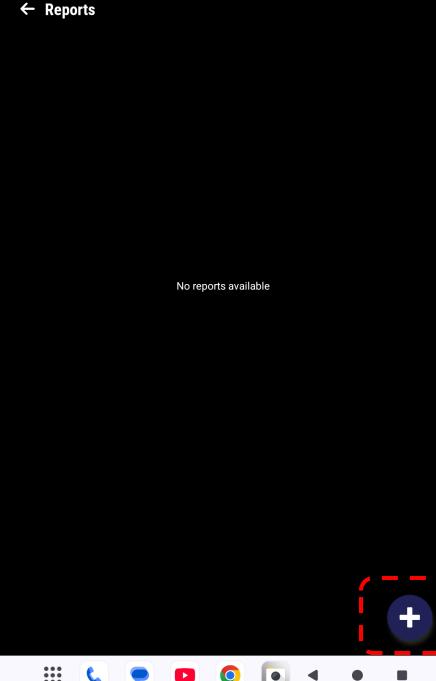






















- The Individual will need to enter their information into the blank spaces and add a logo image
- This information should be remembered by the app for future reports
- Then click "Next"

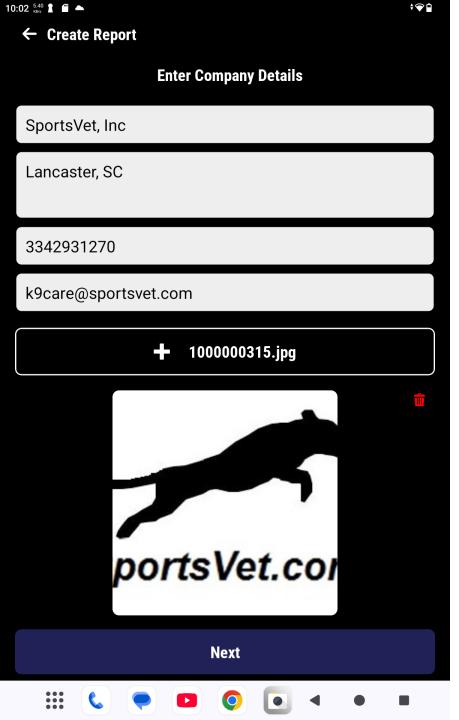


- In this section, the contact information about the subject can be entered
- This section is optional if the user does not want to include this information in the report
- Then click "Next"
- Or the "Back" button if you should want to return to the previous screen

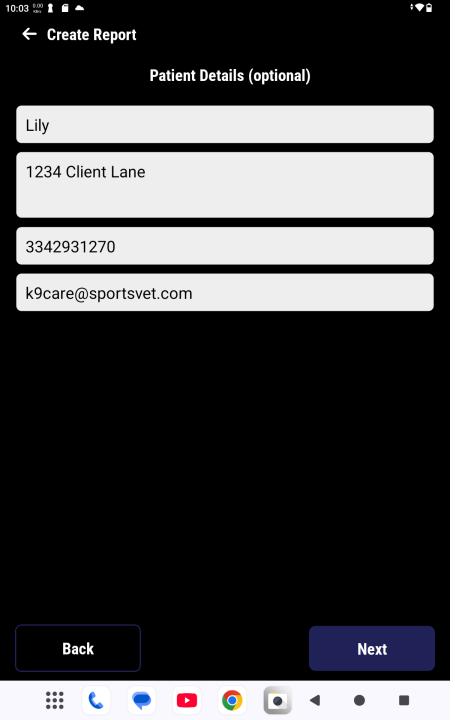

- In this section, the subject's Health information about the subject can be entered
- Date of birth, sex, weight (Kg), Height (cm), Species and Breed
- This section is optional if the user does not want to include this information in the report
- Then click "Next"
- Or the "Back" button if you should want to return to the previous screen

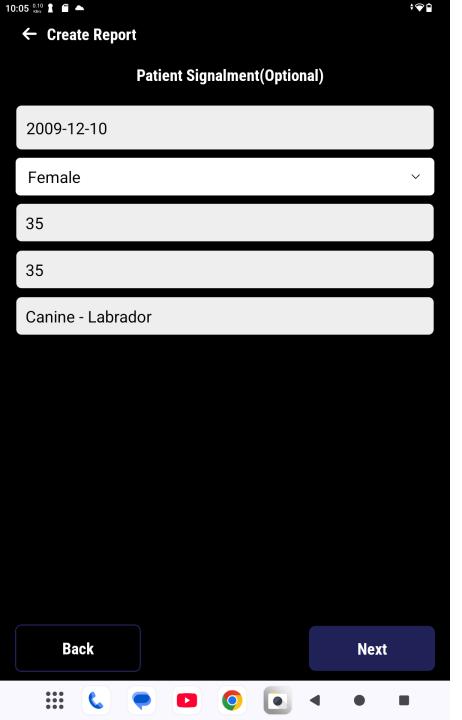


- In this section, general comments about the subject is entered
- This section is "not" optional and will need to be filled in to move to the next section
- Then click "Next"
- Or the "Back" button if you should want to return to the previous screen

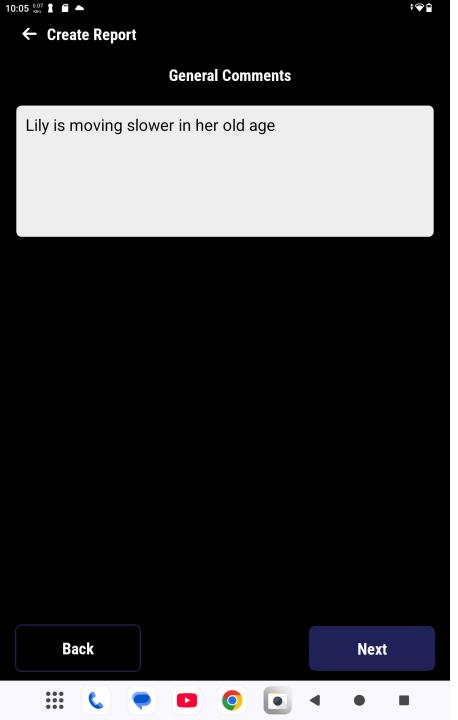





- In this section, a picture is added to the report
- Typically, it is a snapshot with objects and comments that is created in the app
- Comments about the snapshot is entered
- This section is "not" optional and will need to be filled in to move to the next section
- Then click "Next"
- Or the "Back" button if you should want to return to the previous screen

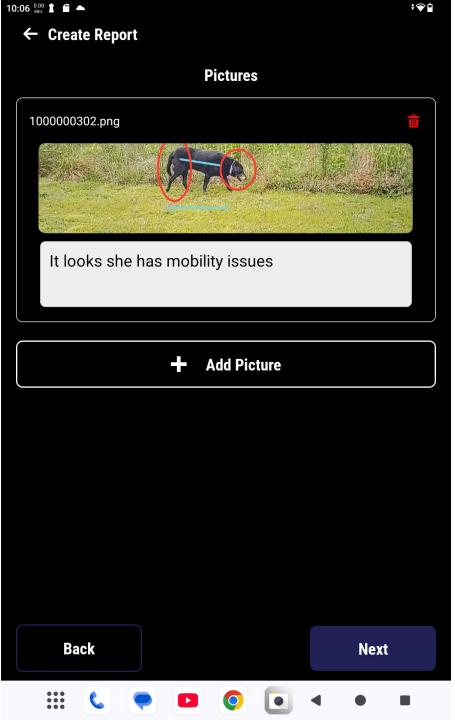

- In this section, other comments are added to the report
- This section is optional if the user does not want to include this information in the report
- Then click "Next"
- Or the "Back" button if you should want to return to the previous screen

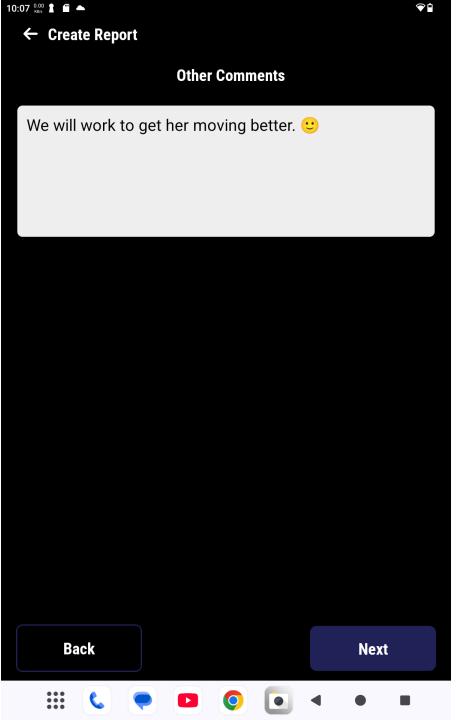

- In this section, specific comments about the analysis and/or plans to manage the subject going forward can be added to the report
- This section is optional if the user does not want to include this information in the report
- Then click "Next"
- Or the "Back" button if you should want to return to the previous screen

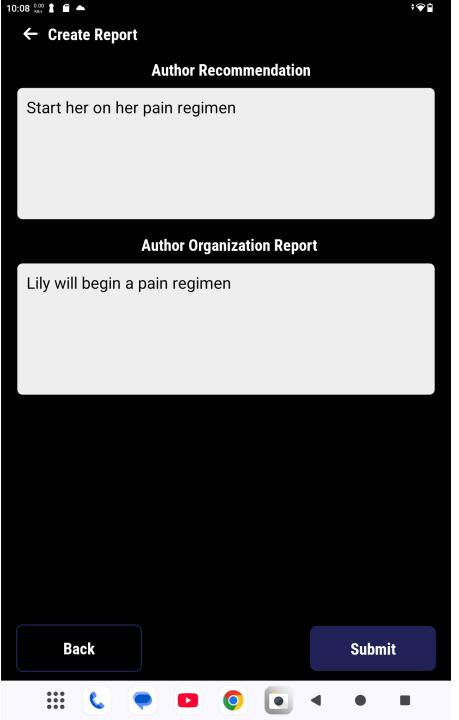











Company Name: SportsVet, Inc

Address: Lancaster, SC

Contact Information: k9care@sportsvet.com | 3342931270


Jul 15, 2025

Patient/Client Name	Client Address	Client Email Address
Lily	1234 Client Lane	k9care@sportsvet.com
Client Mobile Number	Age (DOB)	Sex (M/F)
3342931270	2009-12-10	female
Weight (Kg)	Height (CM)	Species/Breed:
35	35	Canine - Labrador



General Comments

Lily is moving slower in her old age

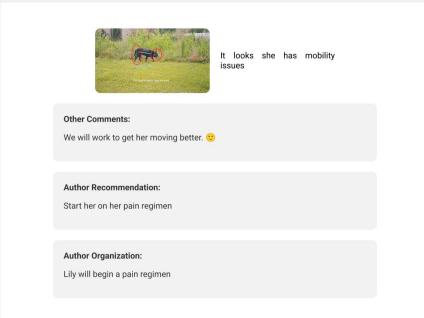

Create Report

- The Report is the generated in a PDF format
- At the top, there are options that allow the user to send a copy of the report to the individual























3/3

- The report will now appear in the "Generated Reports" list
- To create another report, simply click on the "Plus" sign

